



PRESTAN[®]

PRODUCT CATALOG

Prestan® AED Trainers

Prestan® Professional AED Trainer PLUS

The **Prestan Professional AED Trainer PLUS** provides quality and durability like no other trainer in the market. Its size and weight is comparable to an actual AED for a realistic training experience. The AED Trainer is compatible with a remote control.



Patent Pending
Unit Dimensions:
9.5" x 7.25" x 2.75"
Weight is 2.2 lbs.

PRESTAN PROFESSIONAL AED TRAINER PLUS

- PP-AEDT2-101 (Single AED Trainer, English/Spanish)
- PP-AEDT2-102 (Single AED Trainer, English/French)
- PP-AEDT2-401 (4-Pack AED Trainers, English/Spanish)
- PP-AEDT2-402 (4-Pack AED Trainers, English/French)
- PP-APAD2-1 (Single Set of Adult/Child Trainer Pads)
- PP-APAD2-4 (4-Pack Sets of Adult/Child Trainer Pads)

Additional languages available based on distributor-represented geographic area.

Optional Remote Control

PP-AEDT2-100-R

- Simulate loose pad
- Pause, play or stop
- Volume, language and scenario control
- Initiate audible "Press Deeper" during CPR
- Able to control multiple units at once



Prestan® AED UltraTrainer™

The **Prestan AED UltraTrainer** is a compact and customizable unit providing quality, durability and affordability. Sure to be an industry favorite.

PRESTAN AED ULTRATRAINER

- PP-AEDUT-101 (Single AED Trainer, English/Spanish)
- PP-AEDUT-102 (Single AED Trainer, English/French)
- PP-AEDUT-401 (4-Pack AED Trainers, English/Spanish)
- PP-AEDUT-402 (4-Pack AED Trainers, English/French)
- PP-UTPAD-1 (Single Set of dual Adult/Child Pads)
- PP-UTPAD-4 (4-Pack Sets of dual Adult/Child Pads)

Additional languages available based on distributor-represented geographic area.

Patent Pending
Unit Dimensions:
7.25" x 5.25" x 2.0"
Weight is 1.2 lbs.





Prestan® AED TrainerS - features

PRESTAN AED TRAINERS CAN BE CUSTOMIZED FOR ANY TRAINING ENVIRONMENT

The Professional AED Trainer PLUS and the AED UltraTrainer provide the options instructors have been waiting for, combining the most recent technology with realistic features not found on other AED Trainers.

Both AED Trainers are Equipped with Outstanding Features

- Clear, confident voice instructions that provide easy-to-follow directions
- CPR feedback/coaching
- Child defibrillation option
- Training pads are pre-connected and ready to use, unlike other AED trainers that use outdated plug-in technology
- Pre-configured shock scenarios (5 options available)
- Semi-or-fully automatic shocking option
- Compression metronome option (turn on / turn off)
- Rescue Breaths option (turn on / turn off)
- Unit programming can easily toggle between 2 languages
- Programmed with the most up-to-date software and language capability in full compliance with current industry guidelines
- Updates are simple with a module (Professional AED Trainer PLUS) or dongle (AED UltraTrainer)
- Prestan AED Trainers carry a three (3) year warranty

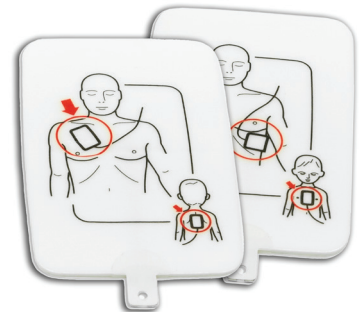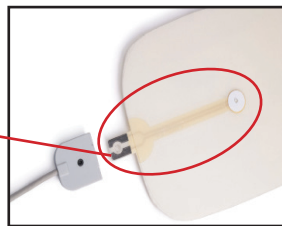

Both AED Trainers are Equipped with Amazing Training Pads

- Pads have an embedded 'switch' which senses proper pad attachment to a manikin. This **Pad Sensing System** (patent pending) allows the Prestan AED Trainers to be used with almost any CPR manikin - no external device needed
- The silicone base adhesive on the pads, coupled with the foam structure, ensures easy application & removal with no residue
- Pads can be used approximately 100 times (with proper care) during CPR/AED classes
- Dual-graphics pads work for both adult and child training

Prestan® Professional Training Manikins

Prestan® Professional Training Manikins

Realistic to the eye and the touch, the **Prestan Professional Training Manikin** series is unlike any other on the market. Available individually or in convenient multi-packs of four for class training purposes, our manikins are uniquely designed to give each student a realistic experience that provides for the best "real life" CPR outcomes.

MANIKIN SINGLES



Prestan Professional Adult Manikin
PP-AM-100M-MS (monitored)
PP-AM-100-MS (non-monitored)



Prestan Professional Child Manikin
PP-CM-100M-MS (monitored)
PP-CM-100-MS (non-monitored)



Prestan Professional Infant Manikin
PP-IM-100M-MS (monitored)
PP-IM-100-MS (non-monitored)

*Manikin single units come with
10 Face-Shield/Lung-Bags
and a nylon Carrying Case*

FEATURES

- Fully compliant with current industry guidelines, as well as the AHA 2019 Integrated Feedback Directive, if purchased (or upgraded) with the CPR Rate Monitor
- Available in Adult, Child and Infant models
- Realistic looking and performing manikins are more acceptable to both instructors and students
- Far easier and faster to carry and setup at training sites
- Patented face/head tilt simulates how an actual victim's head would move if cpr was required and teaches students how to properly open the airway
- Easily inserted face-shield/ lung-bag provides visual chest rise to ensure students receive the most accurate and realistic training
- Face shield is secured so it won't slip and there's no need for adhesives
- Prestan Manikins and parts are not made with natural rubber latex
- Prestan Manikins carry a 3-year warranty

MANIKIN 4-PACKS



Prestan Professional Adult Manikin
PP-AM-400M-MS (monitored)
PP-AM-400-MS (non-monitored)



Prestan Professional Child Manikin
PP-CM-400M-MS (monitored)
PP-CM-400-MS (non-monitored)



Prestan Professional Infant Manikin
PP-IM-400M-MS (monitored)
PP-IM-400-MS (non-monitored)

*Manikin 4-packs come with
50 Face-Shield/Lung-Bags and a nylon Carrying
Case with a detachable Shoulder Strap*



Prestan® Professional Training Manikins

AMAZING INTEGRATED FEEDBACK TOOLS MAKE ALL THE DIFFERENCE!

Prestan's revolutionary *CPR Rate Monitor* provides instant feedback for the rate of chest compressions. Students can gauge speed on their own. Instructors can monitor several students quickly and easily. Lights continually signal rate progression, confirming compressions are within the 100-120 cpm range recommended by the Industry Guidelines.

The yellow light will blink if the 120 cpm upper limit is exceeded, warning students to slow down.



Students can also be confident they are compressing to the appropriate depth. The unique *Clicker Mechanism* provides an audible confirmation when the correct compression depth is reached between 2.0 -2.4 inches.


Clicker is in all Professional Manikins / CPR Rate Monitor is optional



Patented clamshell torso opening allows for quick and easy lung bag insertion and removal

ADULT JAW THRUST MANIKIN

Prestan's Jaw Thrust Head gives instructors the option to teach with the jaw thrust maneuver or head-tilt/chin-lift.

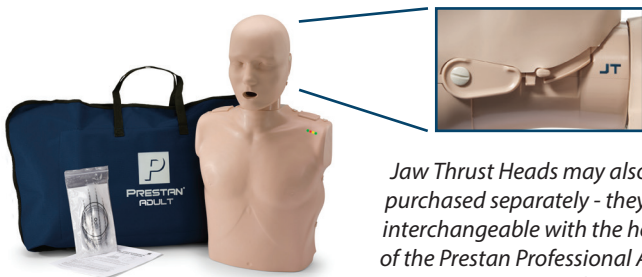

Jaw Thrust Heads may also be purchased separately - they are interchangeable with the heads of the Prestan Professional Adult Training Manikins

Prestan Professional **Adult Jaw Thrust**
Training Manikin

PP-JTM-100M-MS (monitored)
PP-JTM-100-MS (non-monitored)
RPP-JTHEAD-1-MS (head only)

MANIKIN COLLECTIONS

Prestan Professional Manikins with CPR Rate Monitor are available in 2 convenient multi-packs.



The Prestan Collection

(1 Adult / 1 Child / 1 Infant)

PP-FM-300M-MS



The Prestan Family Pack


(2 Adults / 1 Child / 2 Infants)

PP-FM-500M-MS

Multi-packs come with 10 face-shield/lung-bags per manikin, carry case with shoulder strap and pre-installed batteries for the monitors.

UPGRADE YOUR MANIKIN

Easily upgrade existing monitored manikins or install in manikins originally purchased without monitor.

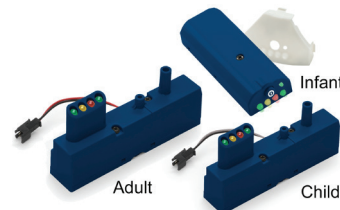

CPR Rate Monitor

for Professional Manikins
CPR Feedback light system includes 120-upper-limit feature (2015 Guidelines)

Adult: RPP-AMON-1
Child: RPP-CMON-1
Infant: RPP-IMON-1



PRESTAN®


Prestan® Ultralite® Manikin 4-Pack

Prestan® Ultralite® Manikins

Prestan is pleased to offer our most portable manikin, the **Ultralite Manikin**. The Ultralite Manikin is amazingly easy to set up and is sold in a convenient and lightweight 4-pack for efficient training-on-the-go. This durable manikin offers an affordable means for CPR training, with all the quality and realism you expect from Prestan.

FEATURES

- Uniquely fast and easy set-up
- Stackable torsos and heads generate a compact package for transporting multiple units
- Clamshell opening provides easy access and lung bag installation
- Manikin is easy to clean
- Compatible with AED Trainer pads
- Realistic head tilt to open airway
- Visible chest rise with lung bag use
- Unique torso web design for realistic feel
- Anatomically correct design includes xiphoid process
- Realistic chest compressions to a 2" depth
- Life tested to a half-million compressions
- Prestan Manikins and parts are not made with natural rubber latex
- Prestan Manikins carry a 3-year warranty

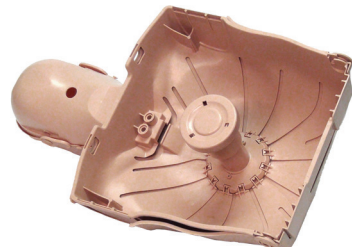

Prestan Ultralite Manikins
PP-ULM-400-MS

Ultralite Manikin 4-pack includes:

- 4 Manikin Torsos
- 4 Manikin Heads
- 4 Compression Pistons
- Nylon Piston Bag
- 50-pack Face-Shield Lung Bags
- Carrying Case
- Detachable Shoulder Strap
- Instruction Sheet



**Complete
4-pack
weighs only
13 pounds!**

Prestan® Products - CONSUMABLES

Prestan® Products - Consumables

DON'T RUN OUT!

Plan ahead to be sure you have enough Prestan Consumables for your training needs. Prestan Consumables are available from your Prestan Authorized Distributor.

FOR YOUR MANIKINS



Face-Shield Lung Bags
Adult 50-Pack: PP-ALB-50
Also available for Child, Infant, and Ultralite Manikins



Rescue Mask Adaptors
10-Bag: 10076-PPA
50-Bag: 10076-PPA-50




Face Shields
Adult 50-Pack: PP-AFS-50
Also available for Child, Infant, and Ultralite Manikins



Generic Face Shield
Adult 36-Count
10077-PPS

FOR YOUR AED TRAINERS

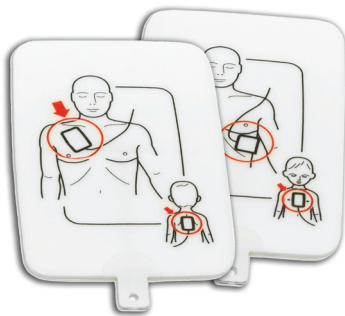

Replacement Training Pads with Adult/Child graphics

For the Professional AED Trainer PLUS:

Single Set: PP-APAD2-1
4-Pack Set: PP-APAD2-4
Replacement Cables: RPP-AEDT2-CABLE

For the AED UltraTrainer:

Single Set: PP-UTPAD-1
4-Pack Set: PP-UTPAD-4
Replacement Cables: RPP-AEDUT-CABLE



Contact your Prestan Authorized Distributor for consumables and replacement parts (not all are listed above).

Consumable items and are not covered by Prestan's three (3) year product warranty.



701 Beta Drive, Suite 3 • Mayfield Village, OH 44143 • 800-434-7989

Visit us at www.prestanproducts.com for more information and product resources

© Prestan Products LLC, January, 2019